




Hipdeer Technology Co.,Ltd

Focused on the manufacturing of medical beauty products

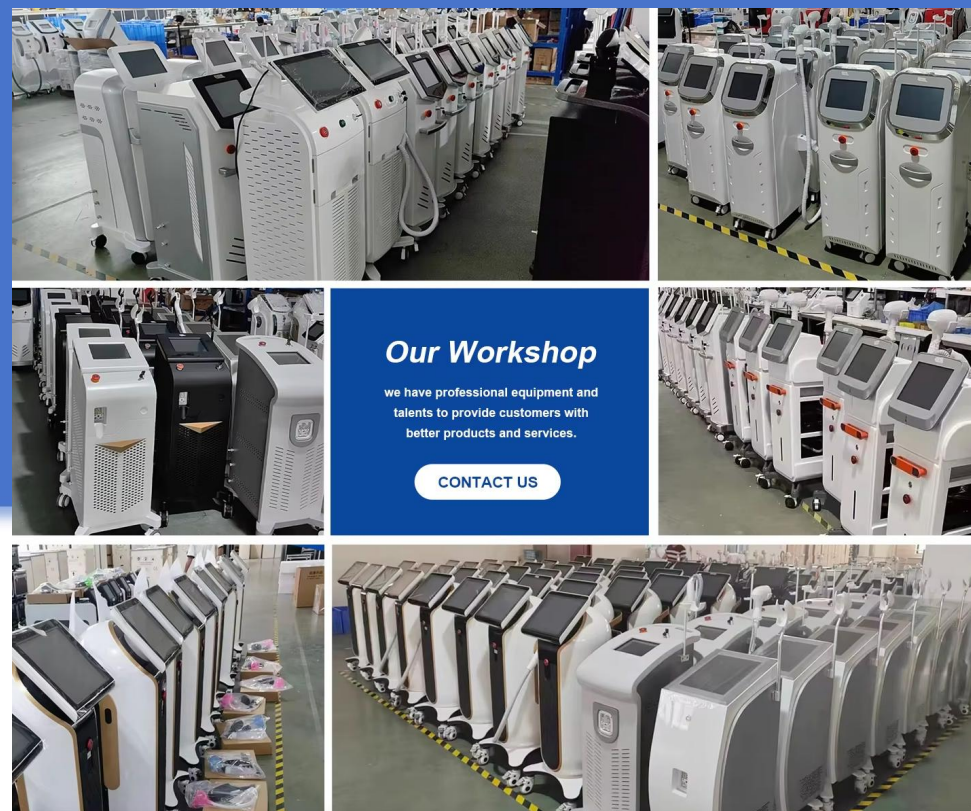

Address: Flat/Rm 745,7/F,Star House,3 Salisbury

Road,Tsim Sha Tsui,Kowloon,Hong Kong

Tel: +86-17815509886

WhatsApp: +86-17815509886

Email: syannsong@hipdeer.com

Online Shop: <https://hipdeer.en.alibaba.com>

PRODUCT CATALOGUE

TECHNOLOGICAL SKINCARE/BEAUTY REGENERATION



Committed to building a world-class enterprise in the beauty device industry and delivering superior products to everyone around the globe.

Technological Innovation & Smart Beauty

PRODUCT CATALOGUE

01/ Magnetic Series

MXS I	01
MXS II	02
Desptop EMS	03
Desptop EMS	04
Mini EMSELLA	05
Private Care Heat Maggie	06
EMface(vertical)	07
EMface(desktop)	08
EMSELLA	09

02/ IPL Hair RemovalSeries

3 in1 IPL	10
2 Screen 3 inl IPL	11
3 in1IPL II	12
Black Whirlwind DPL	13
Pure DPL	14
360 IPL	15
Desktop IPL	16
DPL	17
1.DPL + YAG	18

03/ 808 Hair Removal Series

Lightspeed 808+ Yag	19
Lightspeed 808	20
Sharp 808	21
Sharp 808 + Yag	22
Sharp 808	23
Sharp 808 + Yag	24
Silveredge 808	25
Desktop 808	26
G3000 808+ Yag	27
Korea 808	28
Soprano 808	29
Soprano 808+ Yag	30
Bunny 808	31
Black wirlwind 808	32

04/ Yag Laser Spot Removal Series

Black Whirlwind 808 + Yag	33
Feiton 808	34
Feiton 808 + Yag	35
810 + Yag	36
Mermaid 808	37
Mermaid 808 + Yag	38
Pure 810	39
Pure 810 + Yag	40
Square 808 + Yag	41
Square 808	42
King 808	43
Exclusive 808	44

05/ Pico Laser Series

Pre-ignition Yag	45
Mini Yag	46
Luggage Yag I	47
Luggage Yag II	48
Vertical Yag	49
Desktop Yag	50
Light Speed Yag A	51
Light Speed Yag B	52

Exquisite Pico	53
Picoway	54
Pico Laser Handle	55
Pico Laser Handle	56
Picocare	57
Surperspeed Picosecond	58
Pico Laser Handle	59

06/ Other Series

360 Cryolipolysis	60
CO2 Fractional Pro	61
Eurostar Fractional	62
755	63
Non-invasive Hydrating Glow	64
Exilis Ultra 360	65
Collagen machine	66

MXS I




20,000 muscle contractions in 30 minutes
【Build muscle + burn fat】With waistcoat line, peach butt

02 Handle display



03 Working principle

HIFEM pulses selectively drive motor neurons and it will not affect the skin.
•"Super muscle exercise" stimulates the rapid metabolism of fat cells and powerful decomposition of fat.
•Ultra-fast metabolism causes dysfunction of fat cells and triggers cell apoptosis.
•Has muscle strengthening and growth effects at the same time, "building muscle + burning fat" in one step.

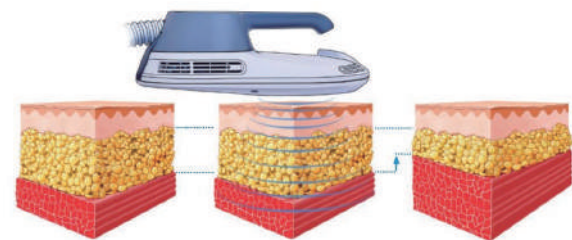

01 Instrument parameters

Total Power Consumption	2800W
GW	68.15kg
Package Size	66*63*153cm
Input Power	220V/110V
Cooling System	Air

MXS II

Technology fat reduction
Painless, Non-invasive

01 Handle display



02 Working principle

HIFEM pulses selectively drive motor neurons and it will not affect the skin.
•"Super muscle exercise" stimulates the rapid metabolism of fat cells and powerful decomposition of fat.
•Ultra-fast metabolism causes dysfunction of fat cells and triggers cell apoptosis.
•Has muscle strengthening and growth effects at the same time, "building muscle + burning fat" in one step.

03 Instrument parameters

Total Power Consumption	5600W
GW	70.3kg
Package Size	71*56*126cm
Input Power	220V/110V
Cooling System	Air



Desptop EMS

4 Handles

NW:38kg
GW:60kg
Power:5600W
Machine Size:58*50*100cm
Package Size:61*64*87cm




2 Handles

NW:24.8kg
GW:42kg
Power:2800W
Machine Size:39*53*45cm
Package Size:57*62*52cm



NW:25.5kg
GW:43kg
Power:2800W
Machine Size:45*45*90cm
Package Size:64*50*57cm


Desptop EMS

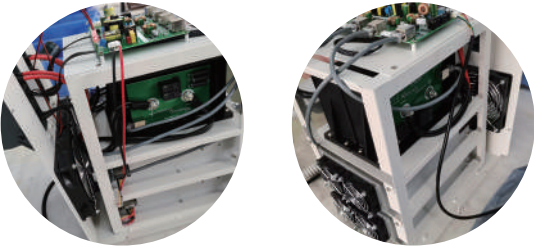

Machine Weight:37.5kg
Power:2800W
Package Size:62*52*52cm



Machine Weight:44.25kg
Power:2800W
Package Size:68.5*53*53cm

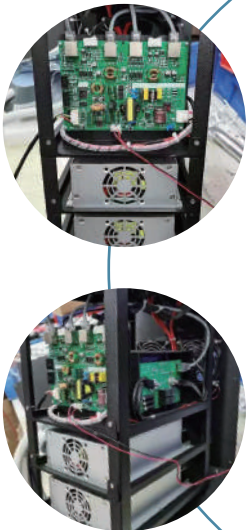

Mini EMSELLA



01 Working principle

Through high-intensity focused electromagnetic fields, motor neurons are stimulated to control the pelvic floor muscles to perform hypercontraction movements, thereby activating the pelvic floor muscle tissue and promoting collagen regeneration. It can effectively exercise the pelvic floor muscles, restore the firmness and elasticity of the pelvic floor muscles, and solve problems such as urinary incontinence, abdominal obesity, and vaginal relaxation

02 Detail display



03 Instrument parameters

Screen Size	4.3 inch
GW	53kg
Package Size	74*74*87 cm
Input Power	220V/110V
Cooling System	Air

Private Care Heat Maggie



02 Working principle

Through high-intensity focused electromagnetic fields, motor neurons are stimulated to control the pelvic floor muscles to perform hypercontraction movements, thereby activating the pelvic floor muscle tissue and promoting collagen regeneration. It can effectively exercise the pelvic floor muscles, restore the firmness and elasticity of the pelvic floor muscles, and solve problems such as urinary incontinence, abdominal obesity, and vaginal relaxation.


01 Detail display



03 Instrument parameters

GW (Machine)	43.35kg
GW(Chair)	63.6kg
Package Size (Machine)	64*63*120cm
Package Size (Chair)	81*77*92cm
Input Power	220V/110V

EMface(vertical)

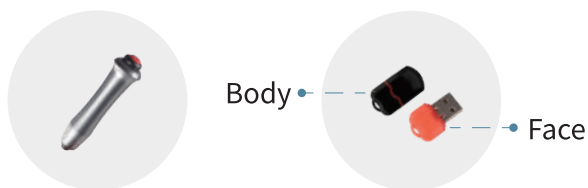

- ① andheld Emergency Stop Rod
- ② USB
- ③ Adapter Box (Use with magnetic tiles)



01 Instrument parameters

Screen Size	12 inch
GW	58.15kg
Package Size	62*48*114cm
Input Power	220V/110V
Cooling System	Air



02 Working principle

MFace working principle: RF and the proprietary HIFES™ technology. HIFES™ is also called high-intensity focused electromagnetic field, which is specially designed to stimulate facial muscles. RF heats skin tissue to stimulate collagen production, thereby improving skin texture and elasticity, The HIFES™ technology achieves tissue remodeling by stimulating deep facial muscles and collagen, effectively improving the aging caused by the loss of facial muscle tissue.

EMface(desktop)



01 Detail display



02 Instrument parameters

Screen Size	12 inch
GW	26.6kg
Package Size	55.5*46*53cm
Input Power	220V/110V
Cooling System	Air



03 Working principle

MFace working principle: RF and the proprietary HIFES™ technology. HIFES™ is also called high-intensity focused electromagnetic field, which is specially designed to stimulate facial muscles. RF heats skin tissue to stimulate collagen production, thereby improving skin texture and elasticity, The HIFES™ technology achieves tissue remodeling by stimulating deep facial muscles and collagen, effectively improving the aging caused by the loss of facial muscle tissue.

EMSELLA



● Orange Pelvic Floor Chair

● Blue Pelvic Floor Chair

01 Working principle

Through high-intensity focused electromagnetic fields, motor neurons are stimulated to control the pelvic floor muscles to perform hypercontraction movements, thereby activating the pelvic floor muscle tissue and promoting collagen regeneration. It can effectively exercise the pelvic floor muscles, restore the firmness and elasticity of the pelvic floor muscles, and solve problems such as urinary incontinence, abdominal obesity, and vaginal relaxation



HIFEMuses 2.5 Tesla intensity and it has no radiation effect on the human body.



Contract and expand muscles 12,000 times in 28 minutes.



Electromagnetic pulse: Stimulates motor neurons to control pelvic floor muscles and causes pelvic floor muscles to contract. It only targets our motor nerves and does not affect other tissues and organs.

02 Detail display



03 Instrument parameters

Screen Size	7inch
GW	63.6kg
Package Size	85.5*74*70 cm
Input Power	220V/110V
Cooling System	Air

3 in1 IPL

01 Handle display



① Yag Laser Handle

Spot removal, Eyebrow tattoo removal, Tattoo removal, Eyeliner tattoo removal, Lip liner tattoo removal, Mole removal, Nevus of Ota removal, Red blood vessels removal, Skin whitening and removing yellowish tone, pore shrinking



② IPL Handle

Painless Ice Point Hair Removal, Whitening Photon Rejuvenation, Pore Contraction, Red Blood Vessel Removal, Acne Removal, Lightening Spots and Freckle Removal



③ RF Handle

Facial wrinkle removal, firming and tightening, smoothing eye-wrinkles, anti-inflammatory, locking in product nutrients.

02 Handle Accessories



IPL Filters



RF Probes



Yag Laser Probes



03 Instrument parameters

Total Power Consumption	2000W
Yag Laser Wavelength	1064/785/532nm
IPL Wavelength	530nm、590nm、640nm
Package Size	61.5*55*133.5cm
GW	43.15kg
Cooling System	Water + Air

2 Screen 3 in1 IPL

01 Handle display



02 Handle Accessories



03 Instrument parameters

Total Power Consumption	2000W
Yag Laser Wavelength	1064/785/532nm
IPL Wavelength	530nm、590nm、640nm
Package Size	63*54.5*130cm
GW	42.55kg
Cooling System	Water + Air

3 in1 IPL II

02 Handle display



03 Handle Accessories

01 Instrument parameters

Total Power Consumption	2000W
Yag Laser Wavelength	1064/785/532nm
IPL Wavelength	530nm、590nm、640nm
Package Size	61.5*55*133.5cm
GW	44.8kg
Cooling System	Water + Air




Black Whirlwind DPL



① DPL Photon Rejuvenation Handle

>Filter Function: 515nm: Used for treating acne and pimples, reducing skin inflammation and scars.560nm: Used for correcting skin pigmentation, reducing dark spots, freckles, and other pigment deposition issues.590nm: Used for enhancing skin radiance, promoting skin regeneration, and reducing fine lines and dullness.615nm: Used for removing spots on the skin, including freckles, sunspots, etc., making the skin more even-toned.640nm: Used for reducing surface redness and vascular dilation on the skin, improving facial redness. 695nm: Used for permanently removing unwanted body hair, achieving smooth skin.

01 filters display

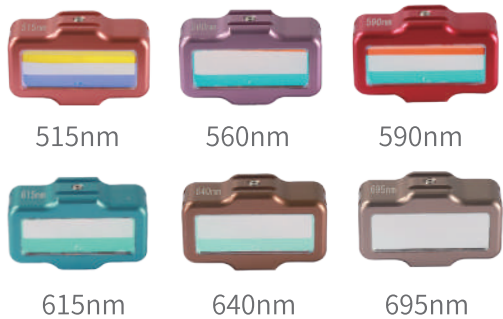

02 Working principle

DPL Precision RejuvenationTechnology is a globally popular high-tech optical cosmetic technology, representing the pinnacle in the field of photon rejuvenation. It utilizes "100-nanometer ultra-narrowband" pulsed light, with targeted energy concentration, enabling rapid and thorough removal of black and red flaws. Simultaneously, it generates photothermal and photochemical effects, rearranging deep collagen and elastin fibers, resulting in smooth and delicate skin and rapid restoration of youthful appearance.

03 Instrument parameters


Total Power Consumption	1200W
GW	60kg
Package Size	61*58.5*122.5cm
Voltage	220V/110V
Cooling System	Water + Air + TEC

Pure DPL



① DPL Photon Rejuvenation Handle

>Filter Function: 515nm: Used for treating acne and pimples, reducing skin inflammation and scars.560nm: Used for correcting skin pigmentation, reducing dark spots, freckles, and other pigment deposition issues.590nm: Used for enhancing skin radiance, promoting skin regeneration, and reducing fine lines and dullness.615nm: Used for removing spots on the skin, including freckles, sunspots, etc., making the skin more even-toned.640nm: Used for reducing surface redness and vascular dilation on the skin, improving facial redness. 695nm: Used for permanently removing unwanted body hair, achieving smooth skin.

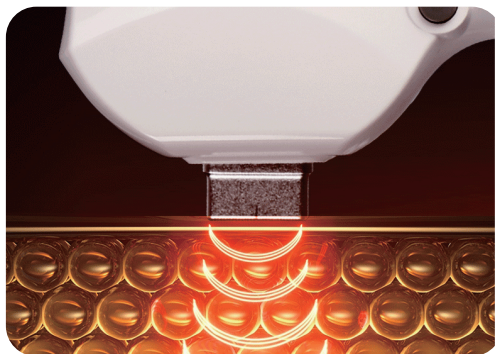

01 Working principle

DPL Precision RejuvenationTechnology is a globally popular high-tech optical cosmetic technology, representing the pinnacle in the field of photon rejuvenation. It utilizes "100-nanometer ultra-narrowband" pulsed light, with targeted energy concentration, enabling rapid and thorough removal of black and red flaws. Simultaneously, it generates photothermal and photochemical effects, rearranging deep collagen and elastin fibers, resulting in smooth and delicate skin and rapid restoration of youthful appearance.

02 Instrument parameters

Total Power Consumption	1200W
GW	60kg
Package Size	61*47*126cm
Voltage	220V/110V
Cooling System	Water + Air + TEC

360 IPL



① Magnetic Light Handle
Painless Ice Point Hair Removal, Whitening Photon Rejuvenation, Pore Contraction, Red Blood Vessel Removal, Acne Removal, Lightening Spots and Freckle Removal

01 Detail display



02 Instrument parameters

Total Power Consumption	1200W
GW	28.7kg
Package Size	60*44*112cm
Voltage	220V/110V
Cooling System	Water + Air




03 Working principle

> Pulse light penetrates the surface layer of the skin, reaching the hair follicle.
> Melanin in the hair follicle and hair shaft absorbs light energy and converts it into heat energy.
> Accumulation of heat energy causes the temperature inside the hair follicle to rise to a certain level (usually exceeding 60 degrees Celsius), leading to the destruction of the hair follicle structure, damage to the hair papilla, and loss of function in the hair stem cells.

Desktop IPL



① Magnetic Light Handle
Painless Ice Point Hair Removal, Whitening Photon Rejuvenation, Pore Contraction, Red Blood Vessel Removal, Acne Removal, Lightening Spots and Freckle Removal

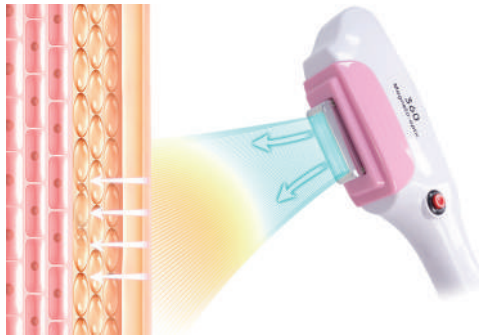

01 Detail display



02 Instrument parameters

Total Power Consumption	1200W
GW	18.75kg
Package Size	59.5*46*41.5cm
Voltage	220V/110V
Cooling System	Water + Air



03 Working principle

> Pulse light penetrates the surface layer of the skin, reaching the hair follicle.
> Melanin in the hair follicle and hair shaft absorbs light energy and converts it into heat energy.
> Accumulation of heat energy causes the temperature inside the hair follicle to rise to a certain level (usually exceeding 60 degrees Celsius), leading to the destruction of the hair follicle structure, damage to the hair papilla, and loss of function in the hair stem cells.

DPL



- ① DPL Photon Rejuvenation Handle
- 2.Filter Function:
480nm/Repair red blood vessels, acne removal, fade acne scars.
530nm/Tender skin, whitening.
590nm/Spot removal.
640nm/Hair removal.
690nm/Hair removal for fine hair.
755nm/Hair removal for dark brown or black skin.



01 Working principle

DPL Precision Rejuvenation Technology is a globally popular high-tech optical cosmetic technology, representing the pinnacle in the field of photon rejuvenation. It utilizes "100-nanometer ultra-narrowband" pulsed light, with targeted energy concentration, enabling rapid and thorough removal of black and red flaws. Simultaneously, it generates photothermal and photochemical effects, rearranging deep collagen and elastin fibers, resulting in smooth and delicate skin and rapid restoration of youthful appearance.

02 Instrument parameters

Total Power Consumption	1200W
GW	37.3kg
Package Size	64*53*121cm
Voltage	220V/110V
Cooling System	Air + Water

DPL + YAG



- ① Yag Laser Handle
- 3.Spot removal, Eyebrow tattoo removal, Tattoo removal, Eyeliner tattoo removal, Lip liner tattoo removal, Mole removal, Nevus of Ota removal, Red blood vessels removal, Skin whitening and removing yellowish tone, pore shrinking



- ② DPL Photon Rejuvenation Handle
- Filter Function:
480nm/Repair red blood vessels, acne removal, fade acne scars.
530nm/Tender skin, whitening.
590nm/Spot removal.
640nm/Hair removal.
690nm/Hair removal for fine hair
755nm/Hair removal for dark brown or black skin



01 Instrument parameters

Total Power Consumption	1200W
GW	40.35kg
Package Size	64*53*121cm
Voltage	220V/110V
Cooling System	Air + Water

02 Working principle

DPL Precision Rejuvenation Technology is a globally popular high-tech optical cosmetic technology, representing the pinnacle in the field of photon rejuvenation. It utilizes "100-nanometer ultra-narrowband" pulsed light, with targeted energy concentration, enabling rapid and thorough removal of black and red flaws. Simultaneously, it generates photothermal and photochemical effects, rearranging deep collagen and elastin fibers, resulting in smooth and delicate skin and rapid restoration of youthful appearance

Lightspeed 808+Pico



① Diode Handle With Touch Screen
Utilizing diode laser with specific wavelength (targeting melanin in hair follicles), for professional hair removal, skin rejuvenation, and whitening, without causing harm to the skin.



② Yag Laser Handle
Spot removal, Eyebrow tattoo removal, Tattoo removal, Eyeliner tattoo removal, Lip liner tattoo removal, Mole removal, Nevus of Ota removal, Red blood vessels removal, Skin whitening and removing yellowish tone, pore shrinking.

② Instrument parameters


Constant Current Power Supply	1500W
GW	35.8kg
Yag Laser Wavelength	532/785/1064nm
Bars	10pcs
Package Size	61.5*56*57cm
Voltage	220V/110V
Language	Support customization

① Working principle
> Using specific wavelengths of 755/808/1064nm achieves better hair removal results with non-irritating emission (operators do not need to wear glasses).
> Shorter hair removal cycles.
> Equipped with a TEC cooling system, the sapphire crystal makes clients comfortable and painless.
> The handle emits light for 50 million pulses, with no cost for consumables

Lightspeed 808



① Diode Handle With Touch Screen
Utilizing diode laser with specific wavelength (targeting melanin in hair follicles), for professional hair removal, skin rejuvenation, and whitening, without causing harm to the skin.

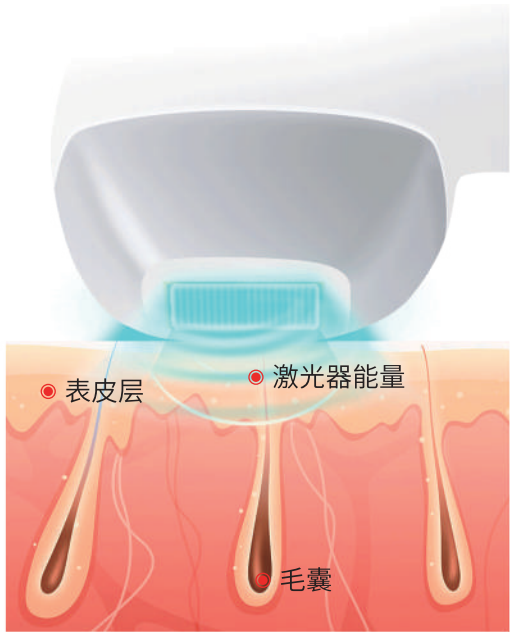

② Instrument parameters


Constant Current Power Supply	1500W
GW	32.35kg
Frequency	1-10HZ (adjustable)
Bars	10pcs
Package Size	61.5*55*57cm
Voltage	220V/110V
Language	Support customization

① Working principle
> Using specific wavelengths of 755/808/1064nm achieves better hair removal results with non-irritating emission (operators do not need to wear glasses).
> Shorter hair removal cycles.
> Equipped with a TEC cooling system, the sapphire crystal makes clients comfortable and painless.
> The handle emits light for 50 million pulses, with no cost for consumables

Sharp 808



① Macro-channel Diode Handle
Utilizing diode laser with specific wavelength (targeting melanin in hair follicles), for professional hair removal, skin rejuvenation, and whitening, without causing harm to the skin.

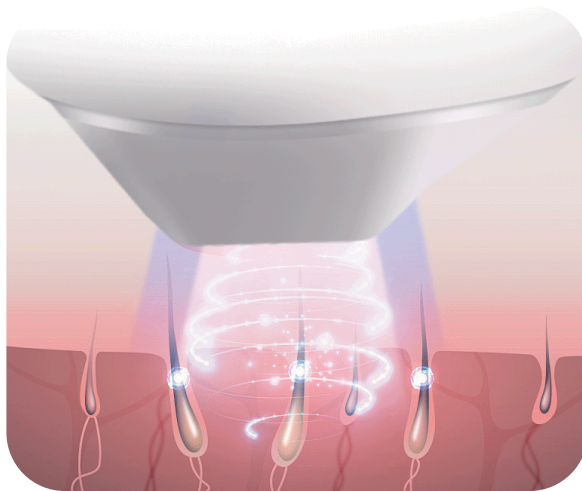


01 Instrument parameters

Bars	6pcs/10pcs
GW	36.85kg
Frequency	1-10HZ (adjustable)
Diode laser wavelength	808nm
Package Size	61.5*55.5*121.5cm
Voltage	220V/110V
Constant Current Power Supply	1500W

02 Working principle

> Using specific wavelengths of 755/808/1064nm achieves better hair removal results with non-irritating emission (operators do not need to wear glasses).
> Shorter hair removal cycles.
> Equipped with a TEC cooling system, the sapphire crystal makes clients comfortable and painless.
> The handle emits light for 50 million pulses, with no cost for consumables

Sharp 808 + Yag

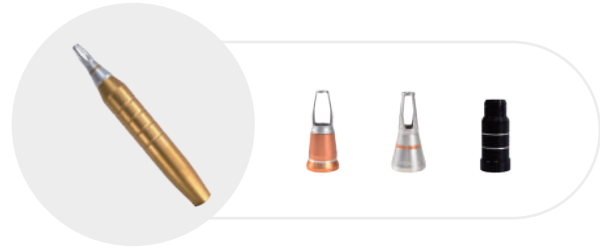

① Yag Laser Handle
Spot removal, Eyebrow tattoo removal, Tattoo removal, Eyeliner tattoo removal, Lip liner tattoo removal, Mole removal, Nevus of Ota removal, Red blood vessels removal, Skin whitening and removing yellowish tone, pore shrinking.



② Macro-channel Diode Handle
Utilizing diode laser with specific wavelength (targeting melanin in hair follicles), for professional hair removal, skin rejuvenation, and whitening, without causing harm to the skin.

01 Internal display



High Quality Wire



Large Water Tank



High Power Supply



24V Diaphragm Pump



02 Instrument parameters


Constant Current Power Supply	1500W
GW	42.65kg
Frequency	1-10HZ (adjustable)
Pulse Width	1-20MS (adjustable)
Package Size	61*55.5*121.5cm
Voltage	220V/110V
Bars	6pcs/10pcs

Sharp 808



① Diode Handle

Pain-free cooling hair removal, skin whitening and rejuvenation, pore refinement, facial redness reduction, freckle and acne removal.

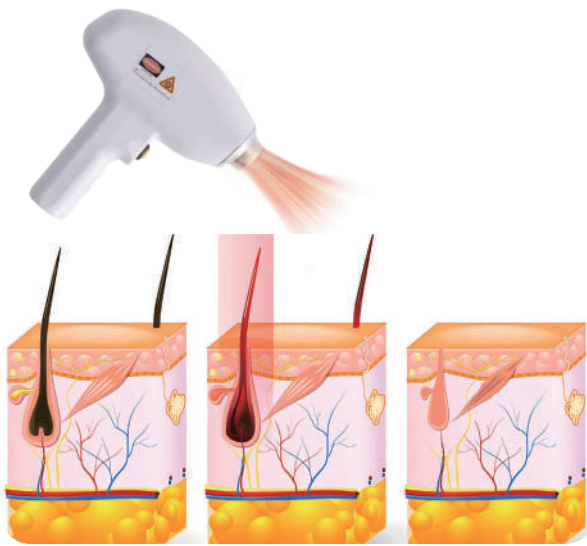

① Instrument parameters


Constant Current Power Supply	1500W
GW	33.9kg
Frequency	1-10HZ (adjustable)
Language	Support Customization
Package Size	61.5*55.5*121.5cm
Voltage	220V/110V
Bars	5pcs



② Working principle

Leveraging the principle of selective photothermolysis, the melanin in hair follicles absorbs light at specific wavelengths (755nm/808nm/1064nm), leading to rapid generation of excessive thermal energy. This targeted heat destroys the hair follicles without damaging the surrounding skin. As the hair follicles are rendered incapable of regrowth, this results in a permanent reduction of hair, achieving the effect of permanent hair removal.

Sharp 808+ Yag

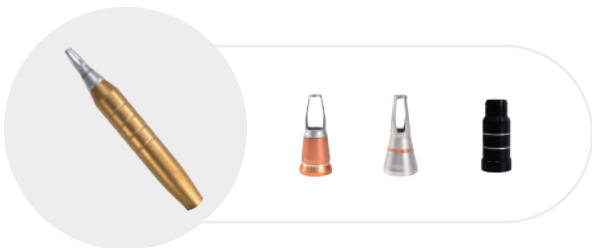

① Yag Laser Handle

Spot removal, Eyebrow tattoo removal, Tattoo removal, Eyeliner tattoo removal, Lip liner tattoo removal, Mole removal, Nevus of Ota removal, Red blood vessels removal, Skin whitening and removing yellowish tone, pore shrinking.



② Diode Handle

Pain-free cooling hair removal, skin whitening and rejuvenation, pore refinement, facial redness reduction, freckle and acne removal.

① Instrument parameters

Constant Current Power Supply	1500W
GW	38.55kg
Yag Laser Wavelength	532/785/1064nm
Language	Support Customization
Package Size	61.5*55.5*121.5cm
Voltage	220V/110V
Bars	5pcs



② Working principle

Leveraging the principle of selective photothermolysis, the melanin in hair follicles absorbs light at specific wavelengths (755nm/808nm/1064nm), leading to rapid generation of excessive thermal energy. This targeted heat destroys the hair follicles without damaging the surrounding skin. As the hair follicles are rendered incapable of regrowth, this results in a permanent reduction of hair, achieving the effect of permanent hair removal.

Silveredge 808



①Macro-channel Diode Handle
Utilizing diode laser with specific wavelength (targeting melanin in hair follicles), for professional hair removal, skin rejuvenation, and whitening, without causing harm to the skin.

01 Detail display



02 Instrument parameters

Constant Current Power Supply	1500W
GW	38.15kg
Frequency	1-10HZ (adjustable)
Language	Support Customization
Package Size	53*49.5*130cm
Voltage	220V/110V
Bars	6pcs/10pcs



03 Working principle

Leveraging the principle of selective photothermolysis, the melanin in hair follicles absorbs light at specific wavelengths (755nm/808nm/1064nm), leading to rapid generation of excessive thermal energy. This targeted heat destroys the hair follicles without damaging the surrounding skin. As the hair follicles are rendered incapable of regrowth, this results in a permanent reduction of hair, achieving the effect of permanent hair removal.

Desktop 808



①Macro-channel Diode Handle
Utilizing diode laser with specific wavelength (targeting melanin in hair follicles), for professional hair removal, skin rejuvenation, and whitening, without causing harm to the skin.

10.4 inch
Touch Screen
Handle Socket
Emergency
Stop Button
Key Switch

Handle Bracket
Energy Conduit
Cart



02 Instrument parameters

Constant Current Power Supply	1500W
GW	28kg
Frequency	1-10HZ (adjustable)
Language	Support Customization
Package Size	61.5*56*57cm
Voltage	220V/110V
Bars	6pcs

01 Working principle

Leveraging the principle of selective photothermolysis, the melanin in hair follicles absorbs light at specific wavelengths (755nm/808nm/1064nm), leading to rapid generation of excessive thermal energy. This targeted heat destroys the hair follicles without damaging the surrounding skin. As the hair follicles are rendered incapable of regrowth, this results in a permanent reduction of hair, achieving the effect of permanent hair removal.

G3000 808+ Yag



① Yag Laser Handle

Spot removal, Eyebrow tattoo removal, Tattoo removal, Eyeliner tattoo removal, Lip liner tattoo removal, Mole removal, Nevus of Ota removal, Red blood vessels removal, Skin whitening and removing yellowish tone, pore shrinking.



② Diode Handle

Pain-free cooling hair removal, skin whitening and rejuvenation, pore refinement, facial redness reduction, freckle and acne removal.

① Instrument parameters

Constant Current Power Supply	1500W
GW	35.2kg
Yag Laser Wavelength	532/785/1064nm
Language	Support Customization
Package Size	58.5*53*113cm
Voltage	220V/110V
Bars	5pcs

② Working principle


Leveraging the principle of selective photothermolysis, the melanin in hair follicles absorbs light at specific wavelengths (755nm/808nm/1064nm), leading to rapid generation of excessive thermal energy. This targeted heat destroys the hair follicles without damaging the surrounding skin. As the hair follicles are rendered incapable of regrowth, this results in a permanent reduction of hair, achieving the effect of permanent hair removal.

Korea 808



① Macro-channel Diode Handle

Utilizing diode laser with specific wavelength (targeting melanin in hair follicles), for professional hair removal, skin rejuvenation, and whitening, without causing harm to the skin.

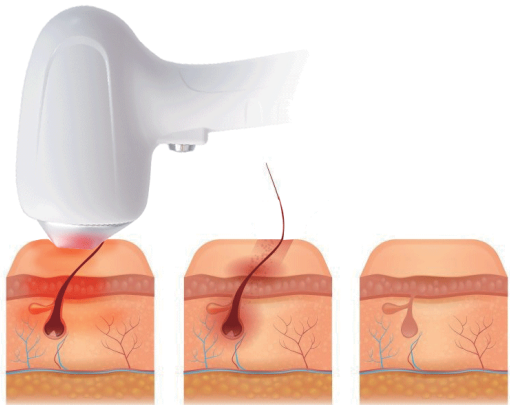

- The wavelength targets and penetrates precisely to the hair follicles.
- Make hair fall out naturally and remove hair permanently
- Stimulate collagen regeneration and shrink pores

① Instrument parameters

Constant Current Power Supply	1500W
GW	42.5kg
Frequency	1-10HZ (adjustable)
Pulse Width	1-20MS (adjustable)
Package Size	65*54*116.5cm
Voltage	220V/110V
Bars	6pcs

② Working principle

> Using specific wavelengths of 755/808/1064nm achieves better hair removal results with non-irritating emission (operators do not need to wear glasses).
> Shorter hair removal cycles.
> Equipped with a TEC cooling system, the sapphire crystal makes clients comfortable and painless.
> The handle emits light for 50 million pulses, with no cost for consumables

Soprano 808



① Diode Handle With Touch Screen
Utilizing diode laser with specific wavelength (targeting melanin in hair follicles), for professional hair removal, skin rejuvenation, and whitening, without causing harm to the skin.

② Internal display



③ Instrument parameters

Constant Current Power Supply	1500W
GW	63.3kg
Frequency	1-10HZ (adjustable)
Pulse Width	1-20MM (adjustable)
Package Size	63*54.5*136cm
Voltage	220V/110V
Bars	12pcs

④ Working principle

Leveraging the principle of selective photothermolysis, the melanin in hair follicles absorbs light at specific wavelengths (755nm/808nm/1064nm), leading to rapid generation of excessive thermal energy. This targeted heat destroys the hair follicles without damaging the surrounding skin. As the hair follicles are rendered incapable of regrowth, this results in a permanent reduction of hair, achieving the effect of permanent hair removal.

Soprano 808+ Yag



① Diode Handle With Touch Screen
Utilizing diode laser with specific wavelength (targeting melanin in hair follicles), for professional hair removal, skin rejuvenation, and whitening, without causing harm to the skin.



② Yag Laser Handle
Spot removal, Eyebrow tattoo removal, Tattoo removal, Eyeliner tattoo removal, Lip liner tattoo removal, Mole removal, Nevus of Ota removal, Red blood vessels removal, Skin whitening and removing yellowish tone, pore shrinking.

③ Instrument parameters

Constant Current Power Supply	1500W
GW	69.6kg
Frequency	1-10HZ (adjustable)
Yag Laser Wavelength	532/785/1064nm
Package Size	63*54.5*136cm
Voltage	220V/110V
Bars	12pcs


④ Working principle

> Using specific wavelengths of 755/808/1064nm achieves better hair removal results with non-irritating emission (operators do not need to wear glasses).
> Shorter hair removal cycles.
> Equipped with a TEC cooling system, the sapphire crystal makes clients comfortable and painless.
> The handle emits light for 50 million pulses, with no cost for consumables

Bunny 808



① Diode Handle With Touch Screen
Utilizing diode laser with specific wave-length (targeting melanin in hair follicles), for professional hair removal, skin rejuvenation, and whitening, without causing harm to the skin.

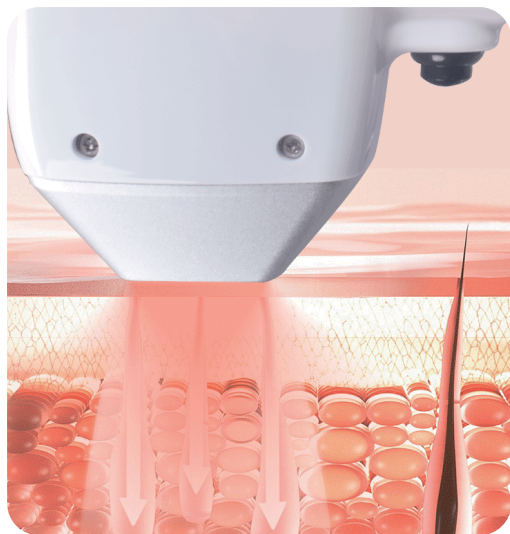


Instrument parameters	
Constant Current Power Supply	1500W
GW	66.7kg
Frequency	1-10HZ (adjustable)
Pulse Width	1-20MM (adjustable)
Package Size	73*48*115cm
Voltage	220V/110V
Bars	12pcs

② Working principle
> Using specific wavelengths of 755/808/1064nm achieves better hair removal results with non-irritating emission (operators do not need to wear glasses).
> Shorter hair removal cycles.
> Equipped with a TEC cooling system, the sapphire crystal makes clients comfortable and painless.
> The handle emits light for 50 million pulses, with no cost for consumables

Black wirlwind 808



① Macro-channel Diode Handle
Utilizing diode laser with specific wave-length (targeting melanin in hair follicles), for professional hair removal, skin rejuvenation, and whitening, without causing harm to the skin

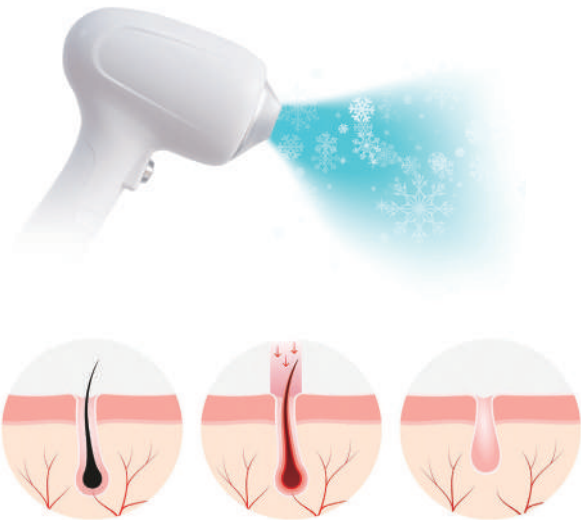


Instrument parameters	
Constant Current Power Supply	1500W
GW	53kg
Frequency	1-10HZ (adjustable)
Pulse Width	1-20MM (adjustable)
Package Size	61*58.5*122.5cm
Voltage	220V/110V
Bars	6pcs/10pcs

② Working principle
> Using specific wavelengths of 755/808/1064nm achieves better hair removal results with non-irritating emission (operators do not need to wear glasses).
> Shorter hair removal cycles.
> Equipped with a TEC cooling system, the sapphire crystal makes clients comfortable and painless.
> The handle emits light for 50 million pulses, with no cost for consumables

Black Whirlwind 808 + Yag



① Macro-channel Diode Handle
Utilizing diode laser with specific wavelength (targeting melanin in hair follicles), for professional hair removal, skin rejuvenation, and whitening, without causing harm to the skin

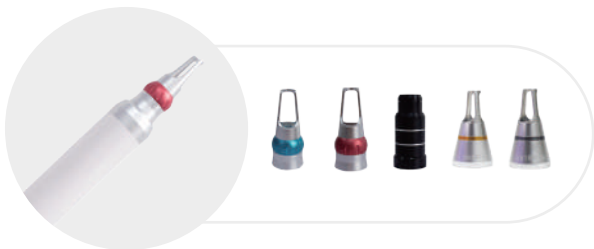

② Yag Laser Handle
5. Spot removal, Eyebrow tattoo removal, Tattoo removal, Eyeliner tattoo removal, Lip liner tattoo removal, Mole removal, Nevus of Ota removal, Red blood vessels removal, Skin whitening and removing yellowish tone, pore shrinking

Instrument parameters

Constant Current Power Supply	1500W
GW	54.1kg
Frequency	1-10HZ (adjustable)
Yag Laser Wavelength	532/785/1064nm
Package Size	61*58.5*122.5cm
Voltage	220V/110V
Bars	6pcs/10pcs

Working principle


>> Using specific wavelengths of 755/808/1064nm achieves better hair removal results with non-irritating emission (operators do not need to wear glasses).
> Shorter hair removal cycles.
> Equipped with a TEC cooling system, the sapphire crystal makes clients comfortable and painless.
> The handle emits light for 50 million pulses, with no cost for consumables.



Feiton 808



① Macro-channel Diode Handle
Utilizing diode laser with specific wavelength (targeting melanin in hair follicles), for professional hair removal, skin rejuvenation, and whitening, without causing harm to the skin

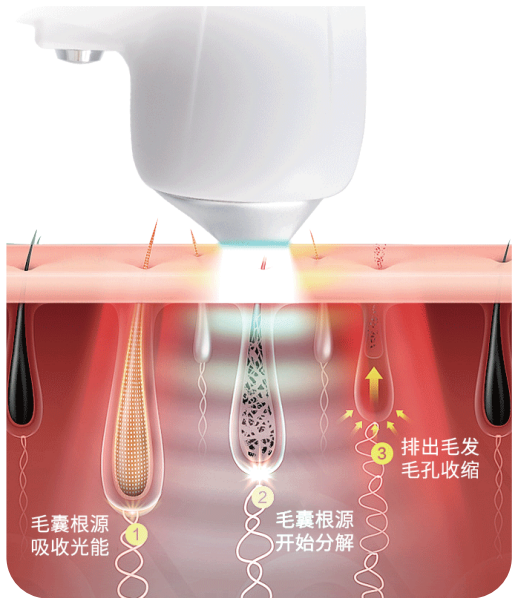

Instrument parameters


Constant Current Power Supply	1500W
GW	51kg
Frequency	1-10HZ (adjustable)
Wavelength	1-20MS (adjustable)
Package Size	71*53.5*121cm
Voltage	220V/110V
Bars	6pcs/10pcs

Working principle

> Using specific wavelengths of 755/808/1064nm achieves better hair removal results with non-irritating emission (operators do not need to wear glasses).
> Shorter hair removal cycles.
> Equipped with a TEC cooling system, the sapphire crystal makes clients comfortable and painless.
> The handle emits light for 50 million pulses, with no cost for consumables.



Feiton 808 + Yag

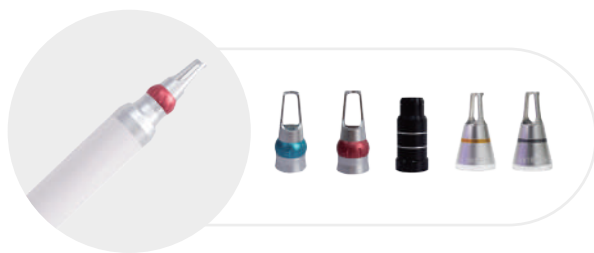

① Yag Laser Handle
Spot removal, Eyebrow tattoo removal, Tattoo removal, Eyeliner tattoo removal, Lip liner tattoo removal, Mole removal, Nevus of Ota removal, Red blood vessels removal, Skin whitening and removing yellowish tone, pore shrinking.



② Macro-channel Diode Handle
Utilizing diode laser with specific wavelength (targeting melanin in hair follicles), for professional hair removal, skin rejuvenation, and whitening, without causing harm to the skin

Instrument parameters

Constant Current Power Supply	1500W
GW	53kg
Yag Laser Wavelength	532/785/1064nm
Frequency	adjustable
Package Size	71*56*126cm
输入电源	220V/110V
巴条	6根/10根



Working principle

> Using specific wavelengths of 755/808/1064nm achieves better hair removal results with non-irritating emission (operators do not need to wear glasses).
> Shorter hair removal cycles.
> Equipped with a TEC cooling system, the sapphire crystal makes clients comfortable and painless.
> The handle emits light for 50 million pulses, with no cost for consumables

810 + Yag



① Macro-channel Diode Handle
3.Utilizing diode laser with specific wavelength (targeting melanin in hair follicles), for professional hair removal, skin rejuvenation, and whitening, without causing harm to the skin.

Detail display



Instrument parameters

Constant Current Power Supply	1500W
GW	52.6kg
Frequency	1-10HZ (adjustable)
Wavelength	1-20MS (adjustable)
Package Size	65*54*116.5cm
Voltage	220V/110V
Bars	6pcs



Working principle

> Using specific wavelengths of 810nm achieves better hair removal results with non-irritating emission (operators do not need to wear glasses).
> Shorter hair removal cycles.
> Equipped with a TEC cooling system, the sapphire crystal makes clients comfortable and painless.
> The handle emits light for 50 million pulses, with no cost for consumables

Mermaid 808



① Macro-channel Diode Handle
Utilizing diode laser with specific wavelength (targeting melanin in hair follicles), for professional hair removal, skin rejuvenation, and whitening, without causing harm to the skin.

01 Detail display



02 Working principle

> Using specific wavelengths of 755/808/1064nm achieves better hair removal results with non-irritating emission (operators do not need to wear glasses).
> Shorter hair removal cycles.
> Equipped with a TEC cooling system, the sapphire crystal makes clients comfortable and painless.
> The handle emits light for 50 million pulses, with no cost for consumables



03 Instrument parameters

Constant Current Power Supply	1500W
GW	37.25kg
Frequency	1-10HZ (adjustable)
Pulse Width	1-20Ms(adjustable)
Package Size	56*51*123cm
Voltage	220V/110V
Bars	5pcs

Mermaid 808 + Yag



① Macro-channel Diode Handle
Utilizing diode laser with specific wavelength (targeting melanin in hair follicles), for professional hair removal, skin rejuvenation, and whitening, without causing harm to the skin.



② Yag Laser Handle
5.Spot removal, Eyebrow tattoo removal, Tattoo removal, Eyeliner tattoo removal, Lip liner tattoo removal, Mole removal, Nevus of Ota removal, Red blood vessels removal, Skin whitening and removing yellowish tone, pore shrinking.



01 Instrument parameters

Constant Current Power Supply	1500W
GW	41.2kg
Frequency	1-10HZ (adjustable)
Yag Laser Wavelength	532/785/1064nm
Package Size	56*50.5*122cm
Voltage	220V/110V
Bars	5 pcs

02 Working principle

> Using specific wavelengths of 755/808/1064nm achieves better hair removal results with non-irritating emission (operators do not need to wear glasses).
> Shorter hair removal cycles.
> Equipped with a TEC cooling system, the sapphire crystal makes clients comfortable and painless.
> The handle emits light for 50 million pulses, with no cost for consumables



Pure 810



① Diode Handle With Touch Screen
Utilizing diode laser with specific wavelength (targeting melanin in hair follicles), for professional hair removal, skin rejuvenation, and whitening, without causing harm to the skin.



Instrument parameters

Constant Current Power Supply	1500W
GW	63kg
Frequency	1-10HZ (adjustable)
Pulse Width	1-20MS (adjustable)
Package Size	61*47*126cm
Voltage	220V/110V
Bars	10pcs

Working principle

> Using specific wavelengths of 810nm achieves better hair removal results with non-irritating emission (operators do not need to wear glasses).
> Shorter hair removal cycles.
> Equipped with a TEC cooling system, the sapphire crystal makes clients comfortable and painless.
> The handle emits light for 50 million pulses, with no cost for consumables

Pure 810 + Yag



① Diode Handle With Touch Screen
Utilizing diode laser with specific wavelength (targeting melanin in hair follicles), for professional hair removal, skin rejuvenation, and whitening, without causing harm to the skin.



② Yag Laser Handle
5.Spot removal, Eyebrow tattoo removal, Tattoo removal, Eyeliner tattoo removal, Lip liner tattoo removal, Mole removal, Nevus of Ota removal, Red blood vessels removal, Skin whitening and removing yellowish tone, pore shrinking.

Instrument parameters

Constant Current Power Supply	1500W
GW	66.1kg
Frequency	1-10HZ (adjustable)
Pulse Width	1-20MS (adjustable)
Package Size	61*47*126cm
Voltage	220V/110V
Bars	10pcs

Working principle

Leveraging the principle of selective photothermolysis, the melanin in hair follicles absorbs light at specific wavelengths (755nm/808nm/1064nm), leading to rapid generation of excessive thermal energy. This targeted heat destroys the hair follicles without damaging the surrounding skin. As the hair follicles are rendered incapable of regrowth, this results in a permanent reduction of hair, achieving the effect of permanent hair removal.

Square 808 + Yag



① Yag Laser Handle
Spot removal, Eyebrow tattoo removal, Tattoo removal, Eyeliner tattoo removal, Lip liner tattoo removal, Mole removal, Nevus of Ota removal, Red blood vessels removal, Skin whitening and removing yellowish tone, pore shrinking.

② Macro-channel Diode Handle
Utilizing diode laser with specific wavelength (targeting melanin in hair follicles), for professional hair removal, skin rejuvenation, and whitening, without causing harm to the skin.

02 Instrument parameters

Constant Current Power Supply	1500W
GW	44kg
Frequency	1-10HZ (adjustable)
Yag Laser Wavelength	532/785/1064nm
Package Size	59*54*122.5cm
Voltage	220V/110V
Bars	5 pcs

01 Working principle

> Using specific wavelengths of 755/808/1064nm achieves better hair removal results with non-irritating emission (operators do not need to wear glasses).

> Shorter hair removal cycles.

> Equipped with a TEC cooling system, the sapphire crystal makes clients comfortable and painless.

> The handle emits light for 50 million pulses, with no cost for consumables.

Square 808



① Macro-channel Diode Handle
Utilizing diode laser with specific wavelength (targeting melanin in hair follicles), for professional hair removal, skin rejuvenation, and whitening, without causing harm to the skin.

01 Detail display



02 Instrument parameters

Constant Current Power Supply	1500W
GW	42kg
Frequency	1-10HZ (adjustable)
Pulse Width	1-20Ms(adjustable)
Package Size	59*54*122.5cm
Voltage	220V/110V
Bars	5 pcs



03 Working principle

> Using specific wavelengths of 755/808/1064nm achieves better hair removal results with non-irritating emission (operators do not need to wear glasses).

> Shorter hair removal cycles.

> Equipped with a TEC cooling system, the sapphire crystal makes clients comfortable and painless.

> The handle emits light for 50 million pulses, with no cost for consumables.


King 808



①



① Macro-channel Diode Handle
Utilizing diode laser with specific wavelength (targeting melanin in hair follicles), for professional hair removal, skin rejuvenation, and whitening, without causing harm to the skin

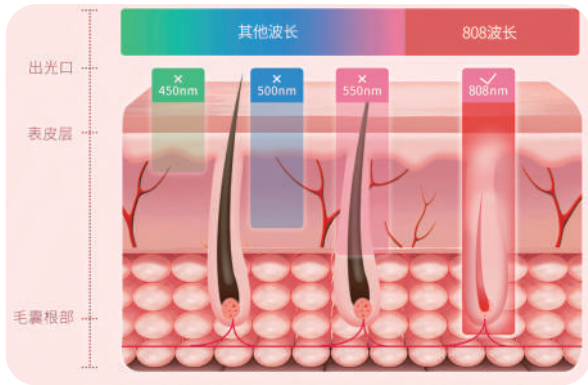


② Instrument parameters

① Working principle

Leveraging the principle of selective photothermolysis, the melanin in hair follicles absorbs light at specific wavelengths (755nm/808nm/1064nm), leading to rapid generation of excessive thermal energy. This targeted heat destroys the hair follicles without damaging the surrounding skin. As the hair follicles are rendered incapable of regrowth, this results in a permanent reduction of hair, achieving the effect of permanent hair removal.

Constant Current Power Supply	1500W
GW	51.2kg
Frequency	1-10HZ (adjustable)
Pulse Width	1-20Ms(adjustable)
Package Size	63.5*53.5*131cm
Voltage	220V/110V
Bars	12 pcs

Exclusive 808

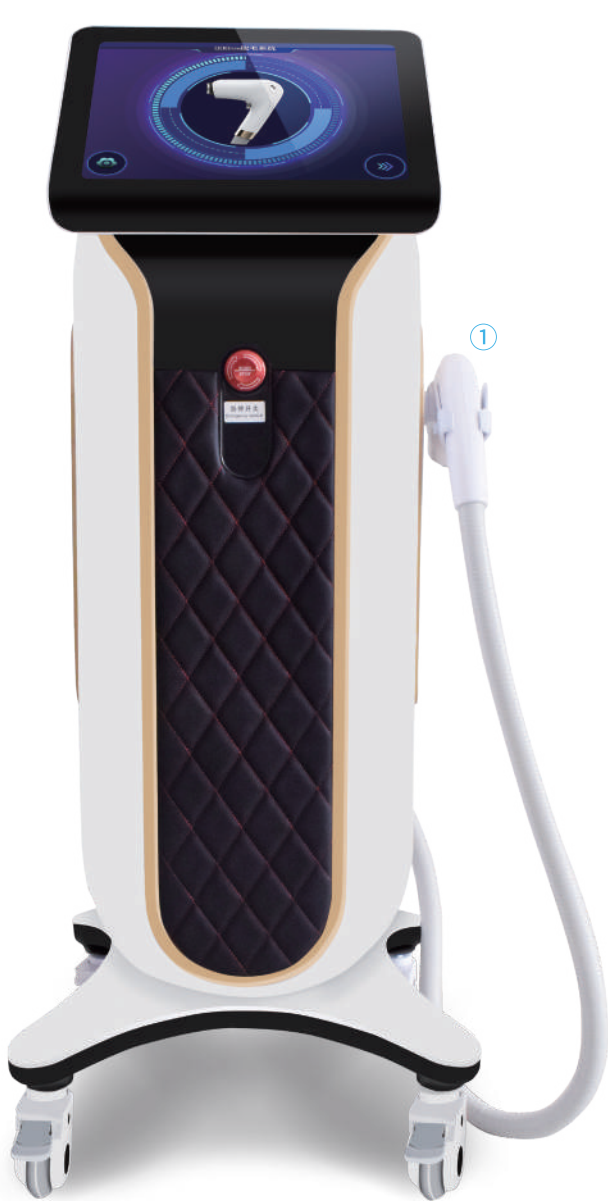


①



① Macro-channel Diode Handle
Utilizing diode laser with specific wavelength (targeting melanin in hair follicles), for professional hair removal, skin rejuvenation, and whitening, without causing harm to the skin

② Internal display

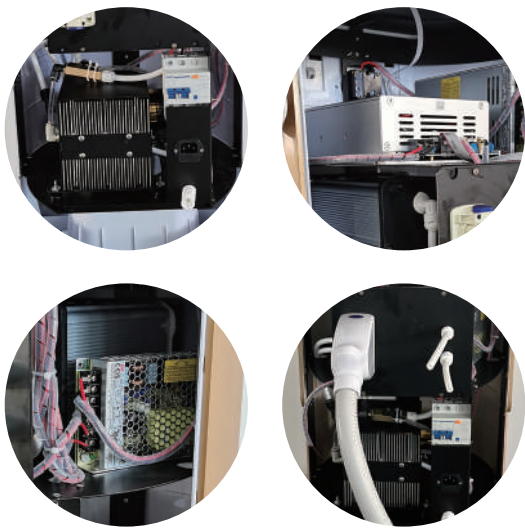


③ Instrument parameters

① Working principle

3.Leveraging the principle of selective photothermolysis, the melanin in hair follicles absorbs light at specific wavelengths (755nm/808nm/1064nm), leading to rapid generation of excessive thermal energy. This targeted heat destroys the hair follicles without damaging the surrounding skin. As the hair follicles are rendered incapable of regrowth, this results in a permanent reduction of hair, achieving the effect of permanent hair removal.

Constant Current Power Supply	1500W
GW	51.85kg
Frequency	1-10HZ (adjustable)
Bars	6 pcs/10 pcs
Package Size	60*55*141.5cm
Voltage	220V/110V
Language	Support Customization

Pre-ignition Yag

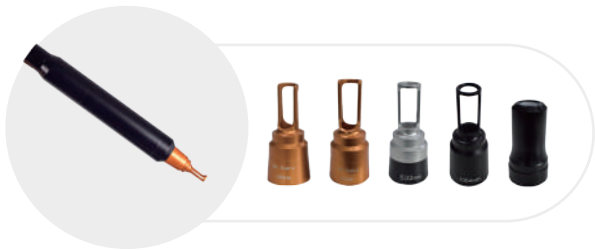

① Yag Laser Handle

Spot removal, Eyebrow tattoo removal, Tattoo removal, Eyeliner tattoo removal, Lip liner tattoo removal, Mole removal, Nevus of Ota removal, Red blood vessels removal, Skin whitening and removing yellowish tone, pore shrinking.

01 Working principle



02 Working principle


The laser emits instantaneously, reaching the chromophores in the diseased area. These chromophores absorb the light energy, heating up rapidly and expanding before instantly shattering into minute particles. Following their ingestion by macrophages, these particles are eventually eliminated from the body via the lymphatic circulation. As a result, the chromophores within the diseased site gradually dwindle until they disappear altogether.



03 Instrument parameters

GW	21kg
Frequency	1-10HZ (adjustable)
Yag Laser Wavelength	532/785/1064nm
Package Size	60*32.5*67cm
Voltage	220V/110V
Language	Support Customization

Mini Yag

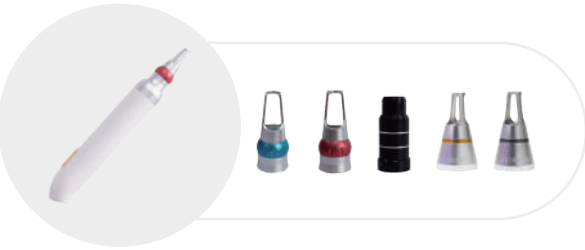

① Yag Laser Handle


Spot removal, Eyebrow tattoo removal, Tattoo removal, Eyeliner tattoo removal, Lip liner tattoo removal, Mole removal, Nevus of Ota removal, Red blood vessels removal, Skin whitening and removing yellowish tone, pore shrinking.

01 Instrument parameters



02 Instrument parameters

Energy	0.5-2.9J/cm ²
GW	16kg
Frequency	1-10HZ (adjustable)
Yag Laser Wavelength	532/785/1064nm
Package Size	59*36*33.5cm
Voltage	220V/110V
Language	Support Customization

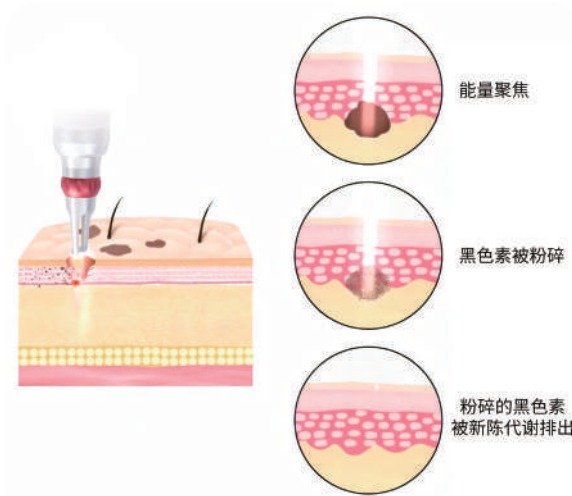

03 Working principle

The laser emits instantaneously, reaching the chromophores in the diseased area. These chromophores absorb the light energy, heating up rapidly and expanding before instantly shattering into minute particles. Following their ingestion by macrophages, these particles are eventually eliminated from the body via the lymphatic circulation. As a result, the chromophores within the diseased site gradually dwindle until they disappear altogether.

Luggage Yag I



① Yag Laser Handle

Spot removal, Eyebrow tattoo removal, Tattoo removal, Eyeliner tattoo removal, Lip liner tattoo removal, Mole removal, Nevus of Ota removal, Red blood vessels removal, Skin whitening and removing yellowish tone, pore shrinking.



01 Instrument parameters



02 Working principle

The laser emits instantaneously, reaching the chromophores in the diseased area. These chromophores absorb the light energy, heating up rapidly and expanding before instantly shattering into minute particles. Following their ingestion by macrophages, these particles are eventually eliminated from the body via the lymphatic circulation. As a result, the chromophores within the diseased site gradually dwindle until they disappear altogether.


03 Instrument parameters

GW	13.1kg
Frequency	1-10HZ (adjustable)
Yag Laser Wavelength	532/785/1064nm
Package Size	60*48.5*45cm
Voltage	220V/110V
Language	Support Customization

Luggage Yag II



02 Instrument parameters

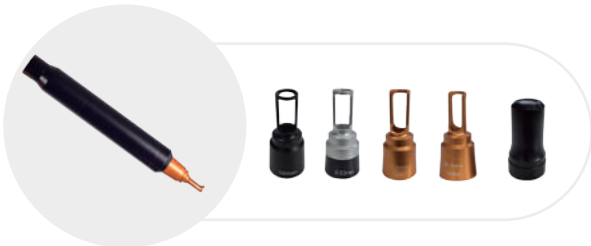

Yag Laser Handle

Spot removal, Eyebrow tattoo removal, Tattoo removal, Eyeliner tattoo removal, Lip liner tattoo removal, Mole removal, Nevus of Ota removal, Red blood vessels removal, Skin whitening and removing yellowish tone, pore shrinking.




01 Instrument parameters

GW	16.5kg
Frequency	1-10HZ (adjustable)
Yag Laser Wavelength	532/785/1064nm
Package Size	54*44*50cm
Voltage	220V/110V
Language	Support Customization

03 Working principle

The laser emits instantaneously, reaching the chromophores in the diseased area. These chromophores absorb the light energy, heating up rapidly and expanding before instantly shattering into minute particles. Following their ingestion by macrophages, these particles are eventually eliminated from the body via the lymphatic circulation. As a result, the chromophores within the diseased site gradually dwindle until they disappear altogether.

Vertical Yag

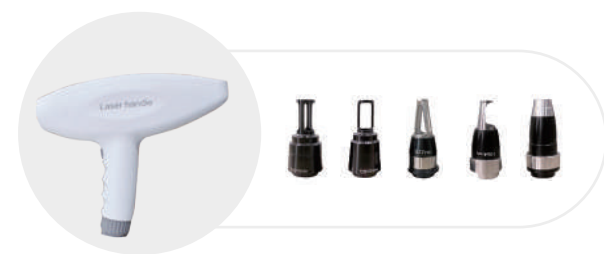

① Yag Laser Handle

Spot removal, Eyebrow tattoo removal, Tattoo removal, Eyeliner tattoo removal, Lip liner tattoo removal, Mole removal, Nevus of Ota removal, Red blood vessels removal, Skin whitening and removing yellowish tone, pore shrinking.

01 细节展示



02 Instrument parameters

GW	36.7kg
Frequency	1-10HZ (adjustable)
Yag Laser Wavelength	532/785/1064nm
Package Size	56*50.5*122cm
Voltage	220V/110V
Language	Support Customization



03 Working principle

The laser emits instantaneously, reaching the chromophores in the diseased area. These chromophores absorb the light energy, heating up rapidly and expanding before instantly shattering into minute particles. Following their ingestion by macrophages, these particles are eventually eliminated from the body via the lymphatic circulation. As a result, the chromophores within the diseased site gradually dwindle until they disappear altogether.

Desktop Yag

① Yag Laser Handle(White Color)

Spot removal, Eyebrow tattoo removal, Tattoo removal, Eyeliner tattoo removal, Lip liner tattoo removal, Mole removal, Nevus of Ota removal, Red blood vessels removal, Skin whitening and removing yellowish tone, pore shrinking.



② Yag Laser Handle(Sliver Color)

Spot removal, Eyebrow tattoo removal, Tattoo removal, Eyeliner tattoo removal, Lip liner tattoo removal, Mole removal, Nevus of Ota removal, Red blood vessels removal, Skin whitening and removing yellowish tone, pore shrinking.



③ Yag Laser Handle(White Color)

Spot removal, Eyebrow tattoo removal, Tattoo removal, Eyeliner tattoo removal, Lip liner tattoo removal, Mole removal, Nevus of Ota removal, Red blood vessels removal, Skin whitening and removing yellowish tone, pore shrinking.



01 Working principle


The laser instantly explodes pigments into tiny particles, making them easier to be metabolized and excreted by the body, thus effectively removing freckles and chloasma.

- 1.Selective photothermal effect
- 2.Pigment decomposition and metabolism
3. Less thermal damage
- 4.Promote collagen regeneration


02 Instrument parameters

Language	Support Customization
White/Black/Red GW	19.1/18.9/18.8kg
Frequency	1-10HZ (adjustable)
Yag Laser Wavelength	562/785/1064nm
Package Size	63.5*44.5*56cm

Light Speed Yag A

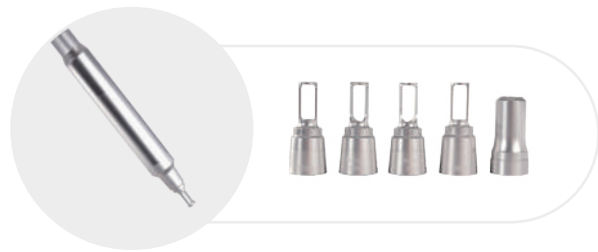

① Yag Laser Handle
Spot removal, Eyebrow tattoo removal, Tattoo removal, Eyeliner tattoo removal, Lip liner tattoo removal, Mole removal, Nevus of Ota removal, Red blood vessels removal, Skin whitening and removing yellowish tone, pore shrinking.

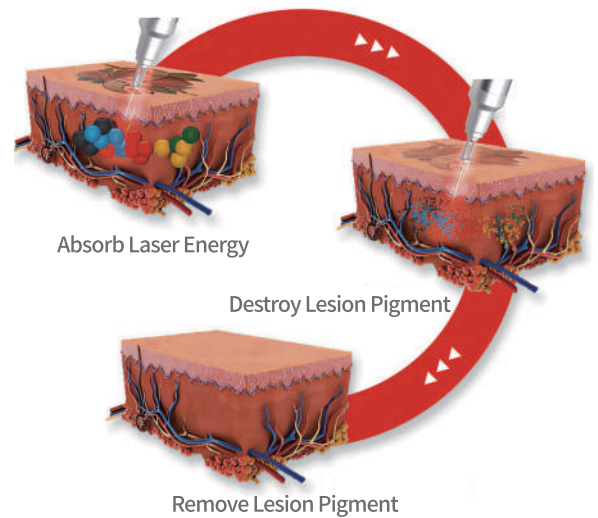

01 Working principle

The laser emits instantaneously, reaching the chromophores in the diseased area. These chromophores absorb the light energy, heating up rapidly and expanding before instantly shattering into minute particles. Following their ingestion by macrophages, these particles are eventually eliminated from the body via the lymphatic circulation. As a result, the chromophores within the diseased site gradually dwindle until they disappear altogether.



02 Instrument parameters

GW	27.95kg
Frequency	1-10HZ (adjustable)
Yag Laser Wavelength	532/785/1064nm
Package Size	61.5*55*57cm
Voltage	220V/110V
Language	Support Customization

Light Speed Yag B



① Yag Laser Handle
Spot removal, Eyebrow tattoo removal, Tattoo removal, Eyeliner tattoo removal, Lip liner tattoo removal, Mole removal, Nevus of Ota removal, Red blood vessels removal, Skin whitening and removing yellowish tone, pore shrinking.

- 532nm
For color pigments
- 无创532nm
Targets color pigments, non-invasive
- 1064nm
Targets dark pigments
- 无创1064nm
Targets dark pigments, non-invasive
- 黑脸娃娃
Skin whitening and rejuvenating

01 Working principle

The laser emits instantaneously, reaching the chromophores in the diseased area. These chromophores absorb the light energy, heating up rapidly and expanding before instantly shattering into minute particles. Following their ingestion by macrophages, these particles are eventually eliminated from the body via the lymphatic circulation. As a result, the chromophores within the diseased site gradually dwindle until they disappear altogether.



02 Instrument parameters


GW	30kg
Frequency	1-10HZ (adjustable)
Yag Laser Wavelength	532/785/1064nm
Package Size	61.5*55*57cm
Voltage	220V/110V
Language	Support Customization

Exquisite Pico



③ Pico Laser Handle

Spot removal, Eyebrow tattoo removal, Tattoo removal, Eyeliner tattoo removal, Lip liner tattoo removal, Mole removal, Nevus of Ota removal, Red blood vessels removal, Skin whitening and removing yellowish tone, pore shrinking.

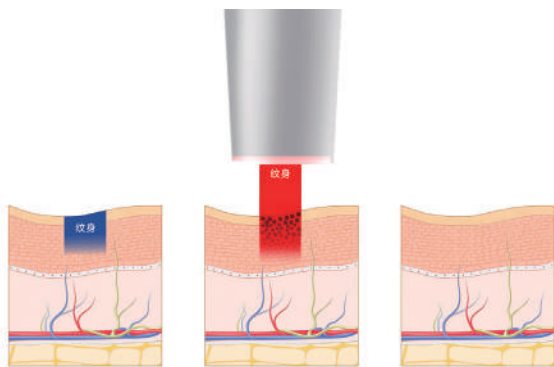

- ①laser arm
- ② hydraulic rod
- ③handle
- ④10.4-inch screen
- ⑤emergency stop button
- ⑥key switch

01 Instrument parameters

Total Power	1500W
GW	92.2kg
Frequency	1-10HZ (adjustable)
Pico Laser Wavelength	532/785/1064nm
Package Size	103*53*135cm
Voltage	220V/110V
Language	Support Customization

02 Working principle


1. Picosecond laser energy is focused 20 times more than traditional lasers, breaking down pigments into fine powder, facilitating the body's metabolism to absorb and expel them.
2. Faster speed allows achieving the desired effect in a very short time, without damaging skin tissues, requiring no recovery period, and providing immediate results.
3. The rapid and powerful energy impact of the picosecond laser causes overwhelming damage to melanin within the skin, preventing the possibility of reactivation.

Picoway



① Pico Laser Handle

Spot removal, Eyebrow tattoo removal, Tattoo removal, Eyeliner tattoo removal, Lip liner tattoo removal, Mole removal, Nevus of Ota removal, Red blood vessels removal, Skin whitening and removing yellowish tone, pore shrinking.

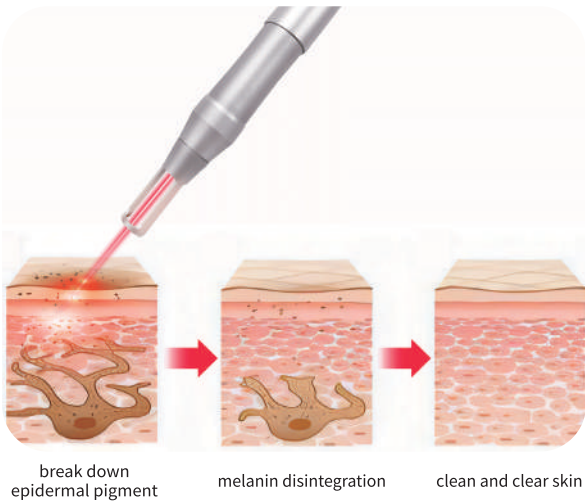

01 Instrument parameters

Total Power	1500W
GW	90.6kg
Frequency	1-10HZ (adjustable)
Pico Laser Wavelength	532/785/1064nm
Package Size	113*51*124.5cm
Voltage	220V/110V
Language	Support Customization

02 Working principle

1. Picosecond laser energy is focused 20 times more than traditional lasers, breaking down pigments into fine powder, facilitating the body's metabolism to absorb and expel them.
2. Faster speed allows achieving the desired effect in a very short time, without damaging skin tissues, requiring no recovery period, and providing immediate results.
3. The rapid and powerful energy impact of the picosecond laser causes overwhelming damage to melanin within the skin, preventing the possibility of reactivation.




Mini super Pico



① Pico Laser Handle

Spot removal, Eyebrow tattoo removal, Tattoo removal, Eyeliner tattoo removal, Lip liner tattoo removal, Mole removal, Nevus of Ota removal, Red blood vessels removal, Skin whitening and removing yellowish tone, pore shrinking.

② DETAIL DISPLAY

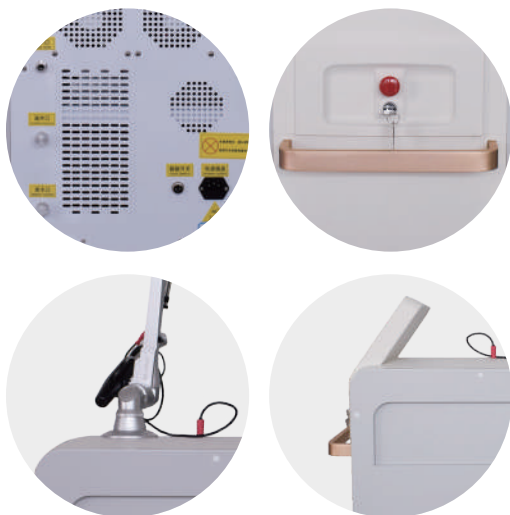

③ Instrument parameters

Total Power	1500W
GW	68.5kg
Frequency	1-10HZ (adjustable)
Pico Laser Wavelength	532/785/1064nm
Package Size	89*46*98cmcm
Voltage	220V/110V
Language	Support Customization

① WORKING PRINCIPLE


1. Picosecond laser energy is focused 20 times more than traditional lasers, breaking down pigments into fine powder, facilitating the body's metabolism to absorb and expel them.
2. Faster speed allows achieving the desired effect in a very short time, without damaging skin tissues, requiring no recovery period, and providing immediate results.
3. The rapid and powerful energy impact of the picosecond laser causes overwhelming damage to melanin within the skin, preventing the possibility of reactivation.

Picosure



① Pico Laser Handle

Spot removal, Eyebrow tattoo removal, Tattoo removal, Eyeliner tattoo removal, Lip liner tattoo removal, Mole removal, Nevus of Ota removal, Red blood vessels removal, Skin whitening and removing yellowish tone, pore shrinking.

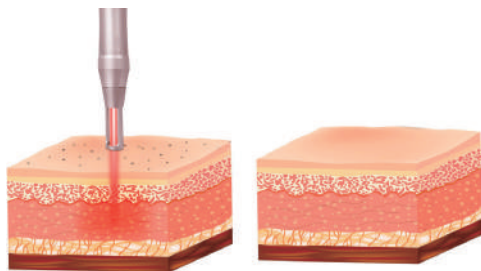

① Instrument parameters



② Instrument parameters

Total Power	1500W
GW	92kg
Frequency	1-10HZ (adjustable)
Pico Laser Wavelength	532/785/1064nm
Package Size	125*61*126cm
Voltage	220V/110V
Language	Support Customization



③ Working principle

1. Picosecond laser energy is focused 20 times more than traditional lasers, breaking down pigments into fine powder, facilitating the body's metabolism to absorb and expel them.
2. Faster speed allows achieving the desired effect in a very short time, without damaging skin tissues, requiring no recovery period, and providing immediate results.
3. The rapid and powerful energy impact of the picosecond laser causes overwhelming damage to melanin within the skin, preventing the possibility of reactivation.

Picosure



① Pico Laser Handle

Spot removal, Eyebrow tattoo removal, Tattoo removal, Eyeliner tattoo removal, Lip liner tattoo removal, Mole removal, Nevus of Ota removal, Red blood vessels removal, Skin whitening and removing yellowish tone, pore shrinking.

② Working principle



② Instrument parameters

Total Power	1500W
GW	94.3kg
Frequency	1-10HZ (adjustable)
Pico Laser Wavelength	532/785/1064nm
Package Size	123*52.5*120cm
Voltage	220V/110V
Language	Support Customization



③ Working principle

1. Picosecond laser energy is focused 20 times more than traditional lasers, breaking down pigments into fine powder, facilitating the body's metabolism to absorb and expel them.
2. Faster speed allows achieving the desired effect in a very short time, without damaging skin tissues, requiring no recovery period, and providing immediate results.
3. The rapid and powerful energy impact of the picosecond laser causes overwhelming damage to melanin within the skin, preventing the possibility of reactivation.

Surperspeed Picosecond



① Pico Laser Handle

Spot removal, Eyebrow tattoo removal, Tattoo removal, Eyeliner tattoo removal, Lip liner tattoo removal, Mole removal, Nevus of Ota removal, Red blood vessels removal, Skin whitening and removing yellowish tone, pore shrinking.

② Working principle



③ DETAIL DISPLAY

Total Power	1500W
GW	76.8kg
Frequency	1-10HZ (adjustable)
Pico Laser Wavelength	532/785/1064nm
Package Size	71*45*115cm
Voltage	220V/110V
Language	Support Customization

② Working principle

1. Picosecond laser energy is focused 20 times more than traditional lasers, breaking down pigments into fine powder, facilitating the body's metabolism to absorb and expel them.
2. Faster speed allows achieving the desired effect in a very short time, without damaging skin tissues, requiring no recovery period, and providing immediate results.
3. The rapid and powerful energy impact of the picosecond laser causes overwhelming damage to melanin within the skin, preventing the possibility of reactivation.



Ultimate Pico



① Pico Laser Handle

Spot removal, Eyebrow tattoo removal, Tattoo removal, Eyeliner tattoo removal, Lip liner tattoo removal, Mole removal, Nevus of Ota removal, Red blood vessels removal, Skin whitening and removing yellowish tone, pore shrinking.

② Detail Display



② Instrument parameters

Total Power	1500W
GW	86.5kg
Frequency	1-10HZ (adjustable)
Pico Laser Wavelength	532/785/1064nm
Package Size	96*47*120cm
Voltage	220V/110V
Language	Support Customization



③ Working principle

1. Picosecond laser energy is focused 20 times more than traditional lasers, breaking down pigments into fine powder, facilitating the body's metabolism to absorb and expel them.
2. Faster speed allows achieving the desired effect in a very short time, without damaging skin tissues, requiring no recovery period, and providing immediate results.
3. The rapid and powerful energy impact of the picosecond laser causes overwhelming damage to melanin within the skin, preventing the possibility of reactivation.

360 Cryolipolysis



② Probe display



② Instrument parameters

Voltage	AC220V±10%、50HZ
GW	70.65kg
Rated Power	800W
RF frequency	5MHz
Package Size	67*66.5*111.5cm
Temperature Range	-10~45°C
Vacuum	0-80Kpa

③ Working principle



Precise and professional negative pressure technology cools the fat cells in the designated area to a specific temperature, causing the triglycerides to change from liquid to solid. Through crystallization and aging, the fat cells gradually die and are excreted from the body through metabolism, achieving the purpose of weight loss. The procedure is quick and recovery is fast, allowing those seeking beauty to achieve more attractive body contours, making them appear well-proportioned and slender.

CO2 Fractional Pro



01 Instrument parameters

GW	71.84kg
Laser Power	60W
Scan Area	20mmx20mm/mj
Package Size	104*53*134.5cm
Scan Speed	>10m/s
Scan Graphic	Circle, triangle, etc. 6 types

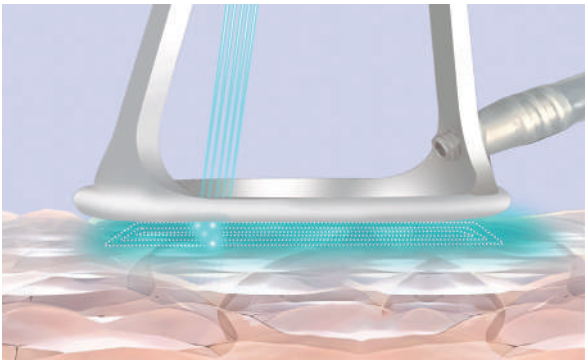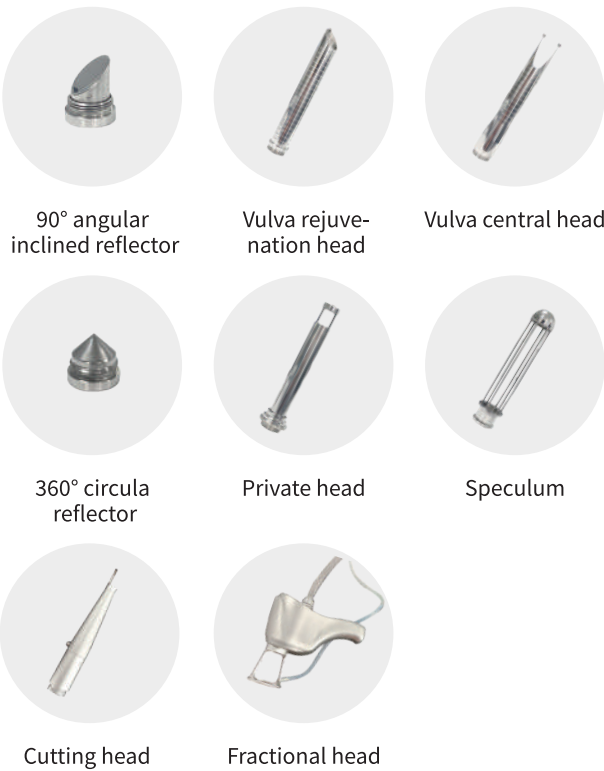

02 Working principle

Fractional laser can emit 50mm-80mm focal spots through a high-focus lens and scan these spots into up to 6 different shapes (circle, square, rectangle, diamond, triangle, line), each suitable for different areas and skin types.

Eurostar Fractional



01 Instrument parameters



02 Instrument parameters

GW	73.15kg
Laser Power	60W
Scan Area	20mmx20mm/mj
Package Size	93*60*121cm
Scan Speed	>10m/s
Scan Graphic	Circle, triangle, etc. 6 types



03 Working principle

Using laser technology, micro-holes with diameters ranging from 100 to 2000 micrometers are evenly created on the skin, triggering a series of biochemical reactions. These micro-holes close within a day, resulting in tightened, rejuvenated skin and wrinkle reduction. Fractional laser operates in a scanning mode with a vibrating mirror. When equipped with a special image generator (CPG), the emission pattern of the light is altered. The fractional laser can emit 50mm-80mm focal spots through a high-focus lens and scan these spots into up to 6 different shapes (circle, square, rectangle, diamond, triangle, line), each suitable for different areas and skin types.




01 Instrument parameters

Total Power	3000W
GW	106.1kg
Spot Size	6-18mm
Wavelength	755/1064nm
Package Size	103*53*135cm
Voltage	220V/110V
Maximum Output Energy	200J/cm ²



① Handle

Using a specific wavelength of 755nm (targeting melanin in hair follicles) for professional hair removal, skin rejuvenation, and whitening, without causing harm to the skin.

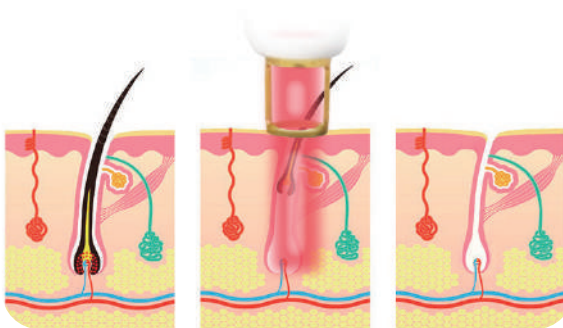

The laser selectively targets the melanin in hair follicles, destroying the hair growth area within the follicle.

Hair naturally falls out, preventing regrowth and achieving the purpose of hair removal.

Hair removal simultaneously stimulates collagen regeneration, reduces pore size, and makes the skin firm and smooth.

02 Working principle


> The long-pulsed 755nm emerald laser meets basic laser hair removal needs. It features an extra-large spot size of 18mm and a repetition rate of 2Hz, resulting in faster and more pronounced results.

> The 755nm wavelength, optimal for melanin absorption, is suitable for Asian individuals with deeper skin tones. It effortlessly removes light-colored and fine hair while whitening the skin and reducing pigmentation.

> The 755nm laser has high melanin absorption. With emerald laser hair removal, there's no need for skin contact or cooling gel, eliminating the risk of cross-infection and skin allergies.

> Dynamic Cooling Device (DCD) ensures treatment safety and comfort, with virtually no pain during the entire procedure.

Non-invasive Hydrating Glow

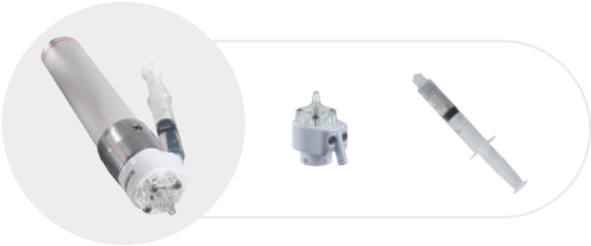

① Non-invasive hydra facial handle


By generating a high-pressure jet stream, the nutrient solution containing extremely small molecules instantly penetrates the epidermis and reaches the basal layer of the skin.

01 Detail display



02 Instrument parameters

Per shot injection volume	Average 0.001ml
Per shot injection frequency	1-20Hz
Per shot injection frequency	0.25mm
Injection depth	100-180v
Nozzle diameter (μm)	100/150/200/250

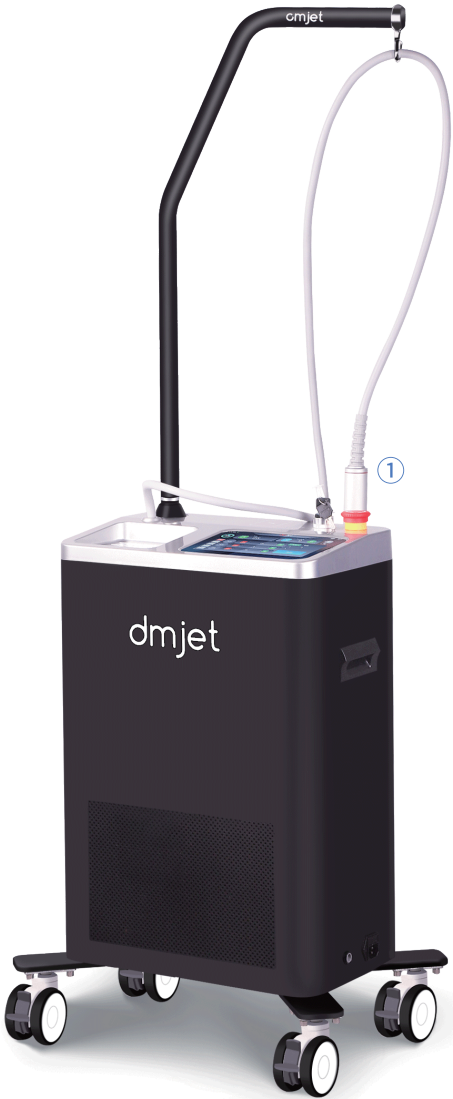

03 Working principle

> High-Pressure Injection Technology: By generating a high-pressure jet stream, this technology instantly penetrates the epidermis with a nutrient solution containing extremely small molecules, delivering it to the basal layer of the skin. This method does not rely on needles to penetrate the skin, avoiding trauma and pain while ensuring the effective absorption of nutrients.

> NanoMicrocrystals: Utilizing high-frequency vibrations or electromagnetic waves, this technology temporarily opens skin channels, allowing the nutrient solution to penetrate more precisely into the target skin layer, enhancing hydration and nourishing effects.

Exilis Ultra 360



① Face Handle
Reduce wrinkles, anti-aging, tighten and lift

② Body Handle
Built-in cooling, fatburning, body shaping

02 Detail display



01 Working principle

Utilizing ultrasound + monopolar radiofrequency technology, its radiofrequency energy precisely reaches different treatment depths within the body, ensuring uniform heat distribution to induce apoptosis of fat cells, thereby achieving fat reduction and body shaping. Its unique collagen regeneration technology effectively stimulates the proliferation of collagen and elastin proteins in the body, reconstructing a new network of collagen fibers to eliminate localized fat accumulation, tighten and lift facial contours, and reduce fine lines. The treatment is safe, reliable, requires no downtime, and prevents rebound, allowing you to achieve skin tightening and body shaping effects while undergoing safe and comfortable treatment, sculpting a perfect figure and reshaping angelic facial features.


03 Instrument parameters

Power Frequency	50HZ
GW	12.2kg
Rated Power	400W
RF Output	150W/40.68HZ
Package Size	48.5*35.5*34.5cm
Voltage	220V/110V
Language	Support Customization

Collagen machine



① RF Handle
Penetrates the dermis layer of the skin up to 5mm, tightens pores, rejuvenates skin, reduces wrinkles, combats aging, and improves complexion.

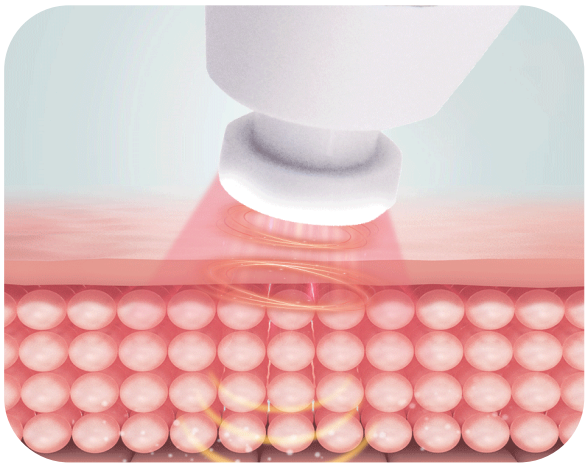

02 Core Advantages

- > Eliminate various wrinkles (nasolabial folds, crow's feet, forehead lines, etc.)
- > Lift marionette lines and enhance facial contours
- > Tighten the outline, reduce baby fat
- > Improve double chin, balance asymmetrical faces
- > Maximum output frequency of 40.68MHz, technologically superior to similar devices
- > Short treatment time, quick results, significant efficacy (immediate effect)
- > Carbon crystal treatment head, corrosion-resistant, high-temperature resistant, no metal crystallization, does not cause skin allergies
- > Adjustable energy levels, temperature adjustable from 38 to 60 degrees
- > Thermal range: subcutaneous 1.5 to 5 cm

01 Instrument parameters

Power Frequency	50HZ
GW	11.25kg
Face Probe	400W
RF Output	150W/40.68HZ Monopolar
Package Size	50*50*31cm
Voltage	220V/110V
Language	Support Customization